



Paint-On Dental Dam

Light-Cure Liquid Dam

Make Lumineers® Clean-up and Finishing, Faster, Easier, Better!




REF 404210

CONTENTS:

- Paint-On Dental Dam syringes, 3 pc. @ 3 ml ea. CE
- Syringe tips (NM-D)
- All-Surface Bonding Chart (NM-D)
- Instruction Card (NM-D)

**Store in cool, dry location
away from light.**

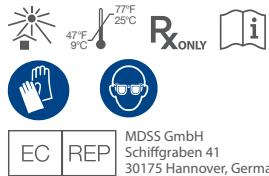

Symbols Glossary available at:
denmat.com/symbols

Den-Mat Holdings, LLC
1017 W. Central Avenue
Lompoc, CA 93436 USA

©2018 Den-Mat Holdings, LLC. All rights reserved.
754402170 R5 11/18CM

Paint-On Dental Dam is a simple, effective and timesaving approach to saliva isolation and soft or hard tissue protection. The low viscosity, colored liquid flows evenly and is easily distinguished from teeth and gingiva. The rubberized resin is easily peeled away when restorations and bleach procedures are completed. Paint-On Dental Dam can also be used for repairing traditional dental dams.

The direct delivery system facilitates restorative procedures and the large capacity (3 ml) syringes cause fewer interruptions.

CAUTION: Wear protective gloves while using this product.

CAUTION: Wear eye protection while using this product.

CAUTION: Paint-On Dental Dam has not been studied in children, pregnant or breast-feeding women.

DIRECTIONS

1. Remove the cap from the syringe (do not discard) and attach the syringe tip to the syringe.
2. Prime the tip by extruding a 1/4-inch line of the material, on a paper pad, before using to ensure adequate dispensing of the material.
3. Apply a bead of Paint-On Dental Dam on gingival surface surrounding the dentition you want to isolate.
4. Carefully manipulate to form a thin film of Paint-On Dental Dam with the syringe tip or spread it with a sable brush.
5. Clean the syringe tip or brush with organic solvent before any ambient light causes polymerization. Then, replace the cap.
6. Hold a light tip 1/2-inch (13 mm) away from the Paint-On Dental Dam material. Light cure, starting with:
 - 5 seconds — Sapphire® Plasma Arc Curing (PAC) Light
 - 10 seconds — Allegro High-Intensity LED Light
 - 30 seconds — Conventional curing light. Repeat this procedure until the entire surface is cured.
7. Peel away the Paint-On Dental Dam using an explorer, Schure 349, or similar instrument after you complete your dental procedure.

STORAGE

Do not expose to temperatures exceeding 77°F (25°C).

Do not expose to direct sunlight.

Do not freeze.

Cold storage is not required.

SAFETY DATA SHEETS AVAILABLE AT denmat.com

TIME-SAVING TIPS

Paint-On Dental Dam is a helpful tool recommended for Lumineers® placement. Before bonding, syringe Paint-On Dental Dam on the lingual side of the interproximals covering the gums and teeth thoroughly. This will make your clean-up easier, faster, better!



Get better impressions with diastemas. Apply Paint-On Dental Dam to interproximals prior to taking the impression. This step will help prevent any tearing of the impression due to the diastema.



Before



After