



ADVANCING 40 YEARS OF CONSISTENT IMPRESSIONS

EXA Advanced™

Hydrophilic Vinyl Polysiloxane Impression Material


EXA Advanced™ Normal & Fast Set:

- Higher tear strength
- Consistently accurate
- No voids, no folds
- Exceptional value, cost effective
- Consistently great impressions every time






Higher Tear Strength Compared to EXAMIX™ & EXAFAST™

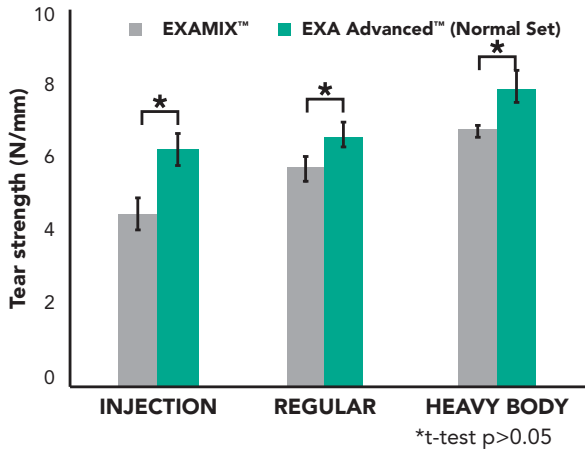


EXA Advanced™ (Normal Set) is an easy-to-use precision VPS impression material that flows better under pressure, is thixotropic, and has working time of up to 2 minutes. It offers exceptional dimensional stability, elastic recovery and tear strength, and “immediate pour” capability.

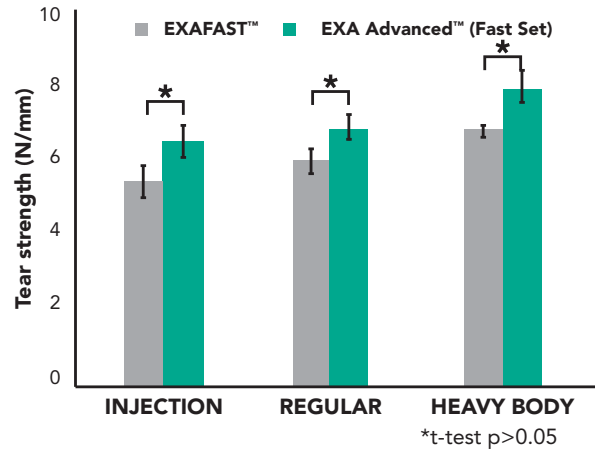

EXA Advanced™ (Fast Set) sets in the mouth in 90 seconds, giving you smooth, precise impressions in a fraction of the normal time for other materials. The short cure time minimizes errors due to patient movement, but provides up to 60 seconds of working time, so you’re not rushed during preparation, making EXA Advanced™ easy to handle and accurate.

Snap Set Cure Rate


Extended Working Time
EXA Advanced™ (Normal Set) has a specifically designed cure rate with a “snap set” feature. This gives you up to 2 minutes to manipulate and seat the material, but with less time in the mouth, something your patients will appreciate.

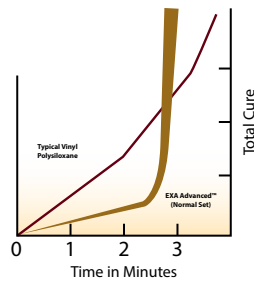

†Data on file at GC America Inc.

Total Cure Time

Faster Impressions = Less Chair Time
With less chair time, you’re more productive and your patients are more comfortable.

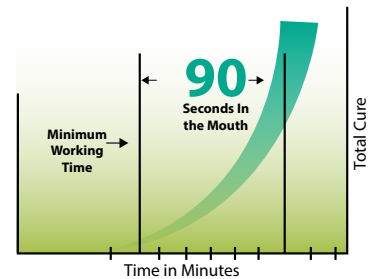

†Data on file at GC America Inc.

Vinyl Polysiloxane Impression Material Usage Guide

Technique	Tray Material	Wash or Syringe Material
Single Mix Also Monophase, Medium Viscosity, Single Phase or Single Impression	Monophase	Monophase
Double Mix	Heavy Body OR Monophase +	Regular OR Injection
Double Arch - Three Function	Monophase OR Heavy Body +	Regular OR Injection
Putty-Wash Also 2-Step Putty Wash	Putty +	Regular Injection
Simultaneous Also 1-Step Putty Wash or Dimple Technique	Putty +	Regular OR Injection

Normal Set:

137116 EXA Advanced™ Heavy Body Normal Set 2 Pack
137117 EXA Advanced™ Regular Type Normal Set 2 Pack
137118 EXA Advanced™ Monophase Normal Set 2 Pack
137119 EXA Advanced™ Injection Type Normal Set 2 Pack
138116 EXA Advanced™ Heavy Body Normal Set 8 Pack
138117 EXA Advanced™ Regular Type Normal Set 8 Pack
138118 EXA Advanced™ Monophase Normal Set 8 Pack
138119 EXA Advanced™ Injection Type Normal Set 8 Pack

Fast Set:

137110 EXA Advanced™ Heavy Body Fast Set 2 Pack
137111 EXA Advanced™ Regular Type Fast Set 2 Pack
137112 EXA Advanced™ Monophase Fast Set 2 Pack
137113 EXA Advanced™ Injection Type Fast Set 2 Pack
138110 EXA Advanced™ Heavy Body Fast Set 8 Pack
138111 EXA Advanced™ Regular Type Fast Set 8 Pack
138112 EXA Advanced™ Monophase Fast Set 8 Pack
138113 EXA Advanced™ Injection Type Fast Set 8 Pack