

Tango-Endo™ Instruments

It Only Takes Two



The Tango-Endo Introductory Kit

Includes: 6 Tango-Endo Instruments (3 Shaping, 3 Finishing), 1-Tango-Endo Reciprocating Handpiece, 60 Precision matched gutta percha points.

See back for details.

Feature

Benefits

Economical:

- Time saving two instrument shaping system.
- Tough and reusable. Use Tango-Endo instruments up to 3 times providing dramatic savings compared to single use systems.
- Affordable Tango-Endo reciprocating handpiece can be used with your current endo motor or any e-type air motor.

Safe:

- Patented flat along the entire length creates two vertical blades, increased cutting, and reduced engagement of the instrument for faster advancement with less resistance to the apex.
- Unique flute design facilitates safe canal shaping and finishing.
- The patented flat increases instrument flexibility without compromising strength.
- Tango-Endo Instruments are specifically designed for the Tango-Endo reciprocating handpiece to virtually eliminate binding and instrument separation.
- Prevents canal distortion.

Efficient:

- Simple 2 step process.
- Superior ability to negotiate curves.
- Quickly navigates to the apex.
- Calibration rings provide quick and precise assistance in determining the working length.



The affordable Tango-Endo reciprocating handpiece can be used with your current endo motor or any e-type air motor.




With Tango-Endo Instruments, it only takes two instruments to complete most cases. By eliminating the fear of fracture you can now complete your endo cases in a simple and predicable process.

Along with an efficient flute design, Tango-Endo instruments have a patented non-interrupted flat along the entire working length.

This design:

- Reduces instrument engagement allowing faster advancement with less resistance.
- Creates space for dentinal debris.
- Reduces compressive and tensile stresses during use enabling improved instrument durability.
- Increased flexibility without sacrificing strength.

The Tango-Endo reciprocation system prevents excessive binding. You can safely work lateral canal walls while preserving canal anatomy.

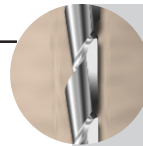

Unique flute design facilitates safe canal shaping and finishing.

ORDERING INFORMATION

Tango-Endo Introductory Kits:

- 5425-00: 6-25mm Tango-Endo Instruments (3 Shaping: 30/02, 3 Finishing:30/04), 1-Tango-Endo Latch-Type Reciprocating Handpiece, 60 Precision matched gutta percha points.
- 5421-00: 6-21mm Tango-Endo Instruments (3 Shaping: 30/02, 3 Finishing:30/04), 1-Tango-Endo Latch-Type Reciprocating Handpiece, 60 Precision matched gutta percha points.

Tango-Endo Instrument Refill Kits:

- 5421-02: Shaping Refill-21mm 6 Tango-Endo Shaping Instruments (30/02)
- 5425-02: Shaping Refill-25mm 6 Tango-Endo Shaping Instruments (30/02)
- 5421-04: Finishing Refill-21mm 6 Tango-Endo Finishing Instruments (30/04)
- 5425-04: Finishing Refill-25mm 6 Tango-Endo Finishing Instruments (30/04)
- 5425-40: Refill-25mm 6 Tango-Endo Instruments for large canals (40/02)
- 5425-50: Refill-25mm 6 Tango-Endo Instruments for extra-large canals (50/02)

Tango-Endo Handpiece

- 5300-00: One Tango-Endo Tango-Endo 4:1 Reduction Latch-Type Reciprocating Handpiece

Tango-Endo Precision Gutta Purcha

- 5400-01: 60 Precision matched gutta percha points, molded (not rolled) to fit your Tango-Endo system.

For our complete line of patented dental products please visit our website at: www.edsdental.com